

Only SomnoMed provides the complete package of products & services for dentists to professionally treat Obstructive Sleep Apnea!



Contact SomnoMed on 1800 445 660 or your local SomnoMed Sales Manager to discuss our trusted Healthcare Partner Referral Network for your patient referrals.

THIS PRODUCT IS NOT AVAILABLE FOR PURCHASE BY THE GENERAL PUBLIC

www.somnomed.com/au



Product Guide



Highest Quality Products, First Time Fit & Stylish Packaging



Evidence: A growing body of published peer reviewed clinical papers.



Guarantee & 3 Year Warranty: Money back guarantee to treat sleep apnea and a comprehensive warranty.



Product Choice: 4 SomnoDent device designs, 2 retention types and additional extras. The SomnoMed product range offers a design suitable for every patient type with a 14 day turn around.



Morning Repositioner: A bite deprogrammer designed to return a patient mandible back to centric occlusion when using a SomnoDent device.



SOMGauge: A superior bite registration tool for accurate & reliable device manufacture.



SomTabs: Specially formulated cleaning tablets.



Referral Networking: SomnoMed consultants will help to link you to local sleep medical providers.



Technical Support: Experienced lab & customer service staff for all technical and product choice enquiries.



Education: Seminars & mentorship programs.



SomnoMed is the **global leader in manufacturing and distributing** oral appliances for the treatment of snoring and sleep apnea.



Building our brand on **heritage, commitment to innovation and performance**, our products are made for comfort and compliance.



Our SomnoDent range of oral appliances has helped over **500,000 patients** worldwide treat their snoring and sleep apnea.

SomnoMed[®] Limited | Level 3, 20 Clarke Street, Crows Nest NSW 2065
Phone: 1800 445 660 | Fax: +61 2 9467 0467 | www.somnomed.com/au
© 2022 SomnoMed Ltd


Clinically Proven
SomnoDent[®]
SLEEP APNEA APPLIANCE



1. Product Range


SomnoDent® SLEEP APNEA APPLIANCE

SomnoDent is the original SomnoMed design consisting of separate upper and lower plates with our patented fin-coupling mechanism. This design is proven to be safe and effective, offering superior comfort with full lip seal and permits normal mouth opening for speech and drinking. SomnoDent is easily adjustable offering 5.0mm forward protrusion and 0.5mm backward from the start position.

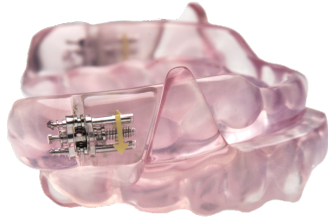

SomnoDent® FUSION™ SLEEP APNEA APPLIANCE

SomnoDent Fusion is an innovation on our original design with the fusion of wing and screw calibration. Titration can be performed precisely with the screw or easily advanced 1mm at a time by changing the wings. This system provides less lateral movement with titration thus more comfort for the patient as their device is advanced. SomnoDent Fusion allows 7.5mm of forward protrusion and 0.5mm backward adjustment, this additional room for titration means less bite resets.

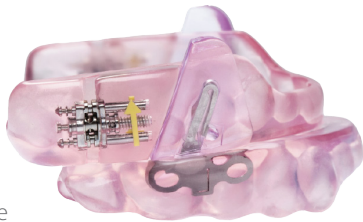

SomnoDent® HERBST ADVANCE™ SLEEP APNEA APPLIANCE

The Herbst Mandibular Advancement Device has long been proven to be efficacious in treating OSA and Snoring. The SomnoMed team have improved the original design by adding a patented visual calibration indicator for objectively measuring protrusion. SomnoDent Herbst Advance can move both laterally and vertically providing superior comfort for heavy grinders where significant lateral freedom is required. Herbst Advance offers up to 7.0mm forward protrusion and 1.0mm backwards from the start position.



2. Retention Type

Flex

Hard Acrylic with Soft Flex Lining

Our unique SMH BFlex provides superior retention and comfort. The Flex material is a soft polymer that is heat cured, packed and pressed; resulting in an intimate fit unlike any other material, anchoring the device even on short teeth with little under cut. With a life expectancy of 3–5 years, Flex is suitable for all dentitions and preferable when crown and bridge work is present.

Classic

Hard Acrylic with Ball Clasp Retention

Made from the highest quality hard acrylic that will not discolour. Classic has been designed to offer durability, comfort and quality. Ball clasps for retention that can be easily adjusted chair side, Classic has life expectancy of 5–7 years.

3. Optional Extras

- Lingual-Less
- Elastic Retention
- EDent (Edentulous)
- Anterior Breathing Space
- Additional Lateral Movement
- DentiTrac Micro-Recorder



Why Choose SomnoDent?

- ✓ Money Back Guarantee*
- ✓ 3 Year Warranty*
- ✓ Easily Adjustable
- ✓ Clinically Validated
- ✓ Custom-fitted
- ✓ Full lip seal
- ✓ Permits normal mouth opening
- ✓ Allows talking and drinking

*Terms & Conditions apply




THE NEW SomnoDent® AVANT™ SLEEP APNEA APPLIANCE

SomnoDent Avant is our first oral appliance that is both digitally designed and manufactured. Both the tough acrylic base and the comfortable BFlex inner liner are milled to create an oral device with unparalleled retention and comfort. The efficacy of the SomnoDent Avant is enhanced by the unique strap mechanism that ensures the lower mandible is extended and the mouth is closed. SomnoDent Avant straps also cater for lateral movement, and exchanging straps allows for forward protrusion of up to 8mm and 1mm backwards from the starting position.



Strap Kit

The strap set is determined by SomnoMed and depend in the patient's jaw size and the bite registration provided by the Dentist. All new SomnoDent AVANT devices are supplied with the designated starting strap assembled on the device and the matching strap kit.

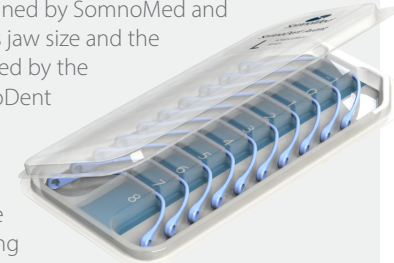

Optional Extras

- Anterior Breathing Hole
A small gap spanning the front teeth, laterally along the occlusal plane may be requested for patient comfort and/or more freedom for forward tongue posturing.
- Anterior Ramp
A small block spanning the two front teeth is built onto the upper tray to disclude the posterior region, i.e. prevent contact at the posterior.
- DentiTrac Micro-Recorder
A micro-recorder may be embedded in the upper splint for the purpose of monitoring patient compliance.